


### OPTIMIZED COMFORT AND SAFETY

EMEDs Soft-Glide® technology has been incorporated in a new family of SCIG needle sets branded under OPTflow®. OPTflow® provides a solution for the increased market demands for faster flow rates and low penetration force to provide optimal comfort and safety to patients and clinicians. With comparable flow rates to 24 gauge SAF-Q® needle sets, OPTflow® needle sets provide optimized comfort. Patented easy-to-use wing design has a stabilizer lip for wing alignment during insertion.

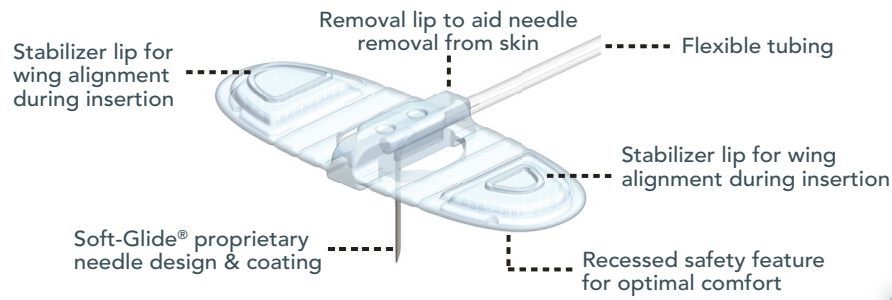

### FEATURES

	SCIG Needle Sets SAF-Q® Models	SCIG Needle Sets OPTflow® Models
• Soft-Glide® coating technology to minimize needle discomfort and anxiety	✓	✓
• The widest selection of single and multi-needle site infusion sets	✓	✓
• Needle lengths: 4, 6, 9, 12 and 14mm to accommodate all age groups and skin types	✓	✓
• Soft translucent wings to facilitate placement and patient comfort	✓	Stabilizer lip for wing alignment to facilitate 90 degree insertion
• Proprietary ultra-flexible polyethylene tubing	✓	✓
• Sterile EMED hypoallergenic site dressing for each needle	✓	Enhanced with dry window
• Proprietary wing design ensures the needle is enclosed completely and securely	✓	Recessed safety feature for optimal comfort
• Gripping lip to aid needle removal from skin	Not a Feature	✓

### BENEFITS

• Needle insertion is virtually painless	✓	Facilitates 90 degree insertion for reduction of pain
• Requires significantly less penetration force when inserting the needle into the skin	✓	✓
• Maximizes patient comfort throughout the infusion	✓	✓
• Great for skin-sensitive patients (thin-skinned, older patients and children)	✓	✓

### ADVANTAGES


• Optimal flow performance with a wide range of fluid viscosities	✓	Maximum flow rate meets and exceeds market requirements
• Easier needle insertion with minimization of tissue damage	✓	✓
• Manufactured with superior polymers to minimize chemical interaction with pharmaceuticals	✓	✓
• Facilitates 90 degree insertion	✓	✓
• Studies show Soft-Glide® administration sets outperform many leading needle manufactures	✓	Optimized to enhance patient experience

\*Comparison report available upon request.


# Increased Safety and Comfort

## With Optimal Flow Performance Aligning with Ever Changing Market Demands


### Quick Reference (See instructions for use for complete guide)

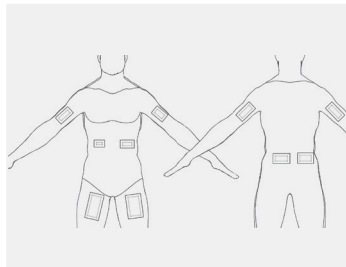


**SELECT** the site of infusion (which should be free from infection). Always choose a site at least 1" to 2" from the previous / other sites.

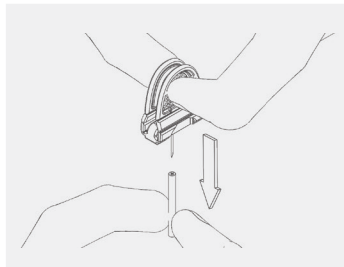

**FOLD** the wings up making sure both wings are aligned and interlocked while holding the wings in between the thumb and index finger and remove needle protector.

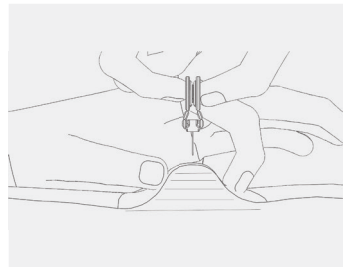

**USE** thumb and index finger to pinch the skin as instructed by your healthcare professional to increase subcutaneous layer thickness. Insert the needle into the skin perpendicular (90°) until the needle is fully inserted into the skin.

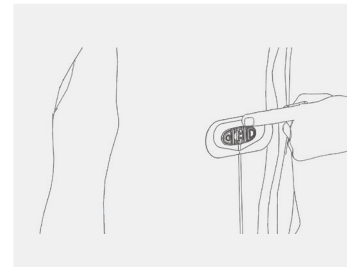

**APPLY** the dressing from over the wings oriented lengthwise and centered to the needle set as close as possible. Then press around the edges of the dressing to fully adhere it to the skin.


### Ordering Information

For orders and inquires, contact your local EMED representative or EMED's customer service department at 916-932-0071. Orders can also be placed via email at [orders@emedtc.com](mailto:orders@emedtc.com).

EMED REORDER NUMBER	TUBING LENGTH	NEEDLE LENGTH (mm)	NEEDLE GAUGE	PRIMING VOLUME	BOX QTY
<b>Single Needle Sets</b>					
OPT12604	70cm	4	26	0.28	50
OPT12606	70cm	6	26	0.28	50
OPT12609	70cm	9	26	0.28	50
OPT12612	70cm	12	26	0.28	50
OPT12614	70cm	14	26	0.28	50
<b>Bi-furcated Needle Sets</b>					
OPT22606	70cm	6	26	0.53	10
OPT22609	70cm	9	26	0.53	10
OPT22612	70cm	12	26	0.53	10
OPT22614	70cm	14	26	0.53	10
<b>Tri-Furcated Needle Sets</b>					
OPT32609	70cm	9	26	0.77	10
OPT32614	70cm	14	26	0.77	10
<b>Quad-Furcated Needle Sets</b>					
OPT42609	70cm	9	26	1.00	10
OPT42612	70cm	12	26	1.00	10
OPT42614	70cm	14	26	1.00	10
<b>Penta-Furcated Needle Sets</b>					
OPT52609	70cm	9	26	1.24	10
<b>Hexa-Furcated Needle Sets</b>					
OPT62609	70cm	9	26	1.47	10

### SOFT-SITE® DRESSING

EMED's patent pending, hypoallergenic dressing is designed for SCIG patient comfort. Window without adhesive avoids unwanted adhesion to needle set.

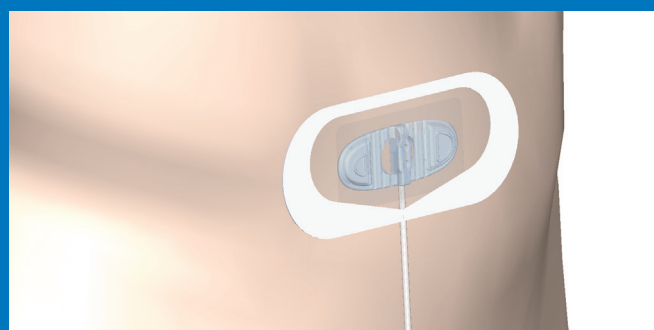

Edges able to be lifted for easy removal to reduce pain and improve patient comfort.

