


Eclipse 5[®]

Designed with your active patients in mind

Key features and benefits

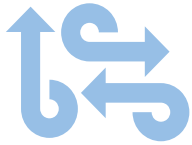

Unmatched Versatility

Setting the standard in transportable oxygen with the broadest oxygen delivery capabilities in continuous and pulse dose. In a compact, transportable design, your patients can go anywhere using the easy-to-attach traveling cart with a telescoping handle and durable wheels. The extended battery life, AC/DC power options and broad portfolio of accessories help to keep your patients active.



Clinically Effective O₂

Designed to adapt to your patients needs with unique clinical features including autoSAT[®], adjustable rise time, UltraSense[™] and myCAIRE[™] to help ensure your patients will get the oxygen they need, when they need it.



Proper O₂ Saturation

With 6 settings up to 3 LPM in continuous flow and 9 settings up to a 192 mL bolus in pulse mode, the Eclipse 5 serves a wide range of patient needs. The unit delivers up to 95% pure oxygen to further ensure proper saturation.



Proven Performance

The Eclipse 5 utilizes Vacuum Pressure Swing Adsorption (VPSA) powerful air separation technology that efficiently eliminates residual nitrogen while delivering the highest purity of oxygen. The robust design features allow for minimum maintenance and long life, providing a device you can rely on.

Specifications

FLOW RATE

Continuous Flow:
6 settings up to 3 LPM

Pulse Flow: 9 settings up
to 192 mL bolus

DIMENSIONS

19.3 x 12.3 x 7.1 in
(49.0 x 31.2 x 18.0 cm)

WEIGHT

15 lbs (without battery)
18.4 lbs (with battery)

SOUND LEVEL

40.0 dB(A) @ Setting 3
Pulse Mode

TRIGGER SENSITIVITY

-0.135-0.37 cm/H₂O
(adjustable)

BATTERY DURATION

up to 5.4 hours in pulse
mode

up to 4.4 hours in
continuous flow mode

POWER

AC Power: 100-240
VAC, 50-60Hz
DC Power: 12 VDC
minimum (all settings
available in AC and DC
power supply)
Battery: Lithium Ion

WARRANTY

3 years
(Includes sieve beds)

1 year on battery and
power supplies



Eclipse 5[®]

Takes you anywhere you want to go

AC POWER SUPPLY (with NEMA power cord)

Item Number
US: 20852326S
EU: 20852325S
UK: 20852324S



UNIVERSAL CART (with telescoping handle)

Item Number
5991-SEQ



ACCESSORY BAG

Item Number
7104-SEQ



DC POWER SUPPLY (with vehicle accessory outlet plug)

Item Number
5942-SEQ



HUMIDIFIER ADAPTER KIT (Humidifier bottle not included)

Item Number
7116-SEQ



TRAVEL CASE

Item Number
5010-SEQ



POWER CARTRIDGE (BATTERY)

Item Number 7082-SEQ



DESKTOP BATTERY CHARGER (Battery Not Included)

Item Number
7112-SEQ



TRAVEL ACCESSORY KIT (Power Cartridge, Desktop Charger, Travel Case)

Item Number
5093-SEQ



PROTECTIVE COVER

Item Number
5052-SEQ



TRAVEL ADAPTER

Item Number 5591-SEQ



WHEELCHAIR PACK

Item Number
5220-SEQ



Find out more about Eclipse 5
and accessory products at

eclipseoxygen.com

SEE PRODUCT WARRANTY STATEMENT FOR COMPLETE INFORMATION.
Please consult the applicable product instructions for use for product indications,
contraindications, warnings, precautions, and detailed safety information.

© 2021 CAIRE Inc. All Rights Reserved. ML-CONC0148 C